IsotopicLabelling R Package: a Practical Guide

Ruggero Ferrazza 2016-03-21

Contents

1	Introduction	1
2	The MS Data Set	1
3	A Compact Way of Processing the Data	2
4	Step-by-Step Processing	3
	4.1 Isotopes and Isotopologues	4
	4.2 Extraction of the Experimental Isotopic Patterns	5
	4.3 Isotopic Pattern Analysis	5
5	Overview of the Results	7
6	Conclusion	8

1 Introduction

The purpose of this document is to explain how to use the *IsotopicLabelling* R package to analyse mass spectrometric isotopic patterns obtained following isotopic labelling experiments.

A typical labelling experiment makes use of substrates enriched in one stable isotope, such as 2 H or 13 C; consequently, after the growth period, there are expected to be metabolites that will have incorporated the labelling isotope, and therefore its relative distribution within those metabolites will be different from its natural occurrence. The IsotopicLabelling R package is based on the principle that, since the isotopic patterns obtained in mass spectrometry reflect the isotopic compositions of the elements making up the observed species, the amount of labelling can be assessed by their proper examination.

Worth of note, because there could be overlapping between the isotopic patterns of different species, the isotopic pattern analysis is better suited for LC-MS or GC-MS data rather than for direct-infusion MS, where the chromatographic step prior to MS detection reduces such issues. Therefore, the current implementation of the package only works for LC-MS or GC-MS data.

Below is a step-by-step explanation on how to use the *IsotopicLabelling* package; the example data set that comes with the package will be used, and the involved functions will be properly introduced and discussed.

2 The MS Data Set

The IsotopicLabelling package requires the MS data to be a data frame with the first two columns representing m/z and retention time (RT) of the mass peaks, and the other columns containing peak intensities or areas (one column for each sample to be analysed).

Since a popular R package for handling MS data is xcms, the IsotopicLabelling R package can also read in xcmsSet objects, which basically contain the peak intensities or areas associated to each of the processed samples, together with their average retention time and m/z. In this case, the function table_xcms is available for converting such objects to the required data frame.

The example data set included in the package is easily accessible:

```
data("xcms_obj")
```

This is an *xcmsSet* object representing lipid extracts of 8 samples from ¹³C labelling experiments:

- The first 4 samples are relative to unlabelled cell cultures (natural ¹³C abundance);
- In the last 4 samples the cells were grown in a substrate where the glucose was replaced by uniformly-labelled ¹³C glucose (99% ¹³C labelling).

This LC-MS data was kindly provided by Dr. Jules Griffin and Dr. Nyasha Munjoma (Department of Biochemistry, University of Cambridge - UK).

The conversion of the *xcmsSet* object to the required data frame is achieved through:

```
peak_table <- table_xcms(xcms_obj)</pre>
```

Here are a few rows of the obtained data frame:

```
rt C12_Sample_1 C12_Sample_2 C12_Sample_3 C12_Sample_4
##
            mz
## 57 157.1582
                                             16463.212
                                                           19478.077
                  39.70566
                              18413.8067
                                                                        18099.8879
## 58 158.1609
                  39.70434
                               1629.4445
                                              1781.991
                                                            1722.311
                                                                         1755.5812
## 59 158.9439 1077.66724
                                599.9926
                                                                          588.8424
                                                    NA
                                                                  NA
      C13_Sample_1 C13_Sample_2 C13_Sample_3 C13_Sample_4
##
                       22425.510
                                     21650.582
## 57
         15949.496
                                                  24001.7999
          1728.988
                        2647.263
                                      1868.273
                                                   1964.7871
## 58
## 59
                 NA
                               NA
                                             NA
                                                    489.4959
```

In addition to *xcms*, this data frame can be obtained in a number of other independent ways, such as through proprietary software of the vendor of the MS instrument; the important point is for the data frame to be properly formatted:

- its first column, named "mz", should contain the mass-to-charge ratios of the peaks;
- its second column, named "rt", should contain the average retention times of the peaks (in seconds);
- the other columns should be named after the samples (one column for sample), and contain peak intensities or areas.

3 A Compact Way of Processing the Data

From the MS data frame, the whole isotopic pattern analysis can be performed through the single, compact function main_labelling. As explained in the reference manual, it requires some input parameters:

• **peak_table**, the data frame containing MS peak intensities or areas;

- compound, a character vector specifying the chemical formula of the compound of interest. A special notation should be used, whereby the character "X" denotes the element with unknown isotopic distribution. For example, the proton adduct of phosphocholine 32:2, [PC 32:2 + H]⁺, has chemical formula C40H77NO8P, but it should be written "X40H77NO8P" for ¹³C labelling experiments, and "C40X76HNO8P" for ²H experiments, in this last case keeping in mind that one hydrogen atom comes from the solvent, and has therefore fixed natural abundance. Please note that adduct ions should be specified, and not the neutral molecular species;
- labelling, a character, either "H" or "C", specifying the labelling isotope;
- mass_shift, the maximum difference between measured and true mass. In other words, the mass accuracy;
- RT, the expected retention time of the compound of interest;
- RT_shift, the maximum difference between true and expected retention time;
- chrom_width, an estimate of the chromatographic width of the peaks;
- initial_abundance, either NA (the default value) or a numeric vector with length equal to the number of samples, containing the initial estimated percentage isotopic abundances of the labelling isotope. If provided, numbers between 0 and 100.

Using the example data set, the parameters to enter for $[PC 32:2 + H]^+$ are:

The output is an object of class *labelling*, a list containing the results of the analysis.

Further details will be given in the following sections, where each of the steps undertaken by main_labelling will be discussed and critically discussed.

4 Step-by-Step Processing

The *IsotopicLabelling* package aims to find the abundance of the isotope used during the labelling experiment within the compounds of interest, based on their MS isotopic patterns. A number of algorithms are available to a-priori compute the isotopic patterns knowing the isotopic abundances; the goal of the *IsotopicLabelling* package is the opposite: starting from measured experimental patterns, this package aims at understanding which is the isotopic distribution giving rise to those patterns. This is basically achieved by a fitting procedure through which the abundance of the labelling isotope is iteratively changed until the best match between predicted and experimental patterns is found.

The main_labelling function basically performs three successive steps:

- 1. It gathers information about the isotopes and all the possible isotopologues arising from the labelling (isotopic_information function);
- 2. It extracts from the experimental data the isotopic patterns of the compound of interest (isotopic_pattern function);
- 3. It analyses the extracted patterns to estimate the percentage isotopic abundance of the labelling isotope.

An alternative to using main_labelling is to run the three distinct functions it is made of, and this route is covered below.

4.1 Isotopes and Isotopologues

[1] "compound" "isotopes" "target"

The first function used by main_labelling is isotopic_information, which summarizes important isotopic information in a single object, a list required by the subsequent functions. The input parameters are the chemical formula and the type of labelling; for $[PC\ 32:2+H]^+$, the list can be obtained through:

```
info <- isotopic_information(compound="X40H77NO8P", labelling="C")</pre>
```

As detailed in the reference manual, the output is a named list:

```
## $names
```

"nTOT"

"nX"

In particular, "isotopes" is a table with the natural isotopic abundances (numbers between 0 and 1) of the elements present in the compound. The two isotopes of the labelling element X are given NA values:

info\$isotopes

```
##
      element
                   mass abundance
## 1
                         0.999885
            Η
               1.007825
## 2
            H 2.014102 0.000115
            N 14.003074 0.996360
## 3
## 4
            N 15.000109
                         0.003640
## 5
            0 15.994915
                         0.997570
            0 16.999132
                         0.000380
## 6
## 7
            0 17.999161
                         0.002050
## 8
            P 30.973762
                         1.000000
## 9
            X 12.000000
                                NA
## 10
            X 13.003355
                                NA
```

Importantly, "target" is a named vector with the exact masses of all the possible isotopologues arising because of the labelling isotope; in the example, [PC 32:2 + H]⁺ has 40 carbon atoms, and therefore the possible isotopologues coming from 13 C span a 41 mass range: the lightest one is the monoisotopic species (with 40 12 C atoms), whereas the heaviest is the species with 40 13 C atoms. However, the isotopic patterns also depend on the other elements, and therefore the list of target isotopologues is further extended by two m/z units, enough for small and medium-sized molecules such as lipids and metabolites.

The naming of the target masses follows this convention:

- M+0 is the monoisotopic mass, the sum of the masses of the atoms using the lightest isotope for each element, X included;
- M+1 is the mass where one light isotope (either X or any other element) is replaced by its heaviest counterpart;
- M+i is the mass where there have been "i" replacements.

The underlying assumption is that the MS resolution is not high enough to resolve the isotopic fine structure; consequently, the replacement of, for example, ¹H with ²H is indistinguishable from ¹²C with ¹³C. This is true for most of the instruments currently used in LC-MS measurements.

4.2 Extraction of the Experimental Isotopic Patterns

Once the target masses are known, the experimental isotopic patterns are extracted from the MS data through the function isotopic_pattern. Keeping on with the example of [PC 32:2 + H]⁺, this can be achieved through the command:

In addition to the table of peaks and the list obtained above, the used also has to provide the parameters mass_shift, RT, RT_shift and chrom_width (already discussed above and further detailed in the reference manual).

The output is a matrix containing the extracted signals, with its first two columns reserved for the exact m/z and the retention times of the peaks:

```
head(experimental_patterns, n=3)
```

```
##
                       rt C12_Sample_1 C12_Sample_2 C12_Sample_3 C12_Sample_4
             mz.
## M+0 730.5387 282.8571
                              2033683.5
                                           1504910.9
                                                         1702717.5
                                                                       1556318.5
## M+1 731.5420 282.8571
                               918440.7
                                             682545.9
                                                          772604.1
                                                                        706638.7
## M+2 732.5454 282.8571
                               219390.3
                                             169205.2
                                                          189325.5
                                                                        169515.6
       C13_Sample_1 C13_Sample_2 C13_Sample_3 C13_Sample_4
##
## M+O
                   0
                                 0
                                               0
                                                            0
                   0
                                 0
                                               0
                                                            0
## M+1
                                               0
## M+2
                   0
                                 0
                                                             0
```

Each of its columns, therefore, represents the extracted experimental pattern for that sample. To get this result, the isotopic_pattern function does the following:

- 1. For all the masses in the "target" vector, it finds and stores the indices of the peaks in the experimental data frame that are within the specified m/z and RT ranges;
- 2. It considers the retention times of the obtained indices, and compares them across isotopologues in order to group them: peaks that differ in RT within the specified chromatographic width are assumed to be two isotopologues;
- 3. More groups may have been identified: those containing less than two isotopologues are discarded and, if still more groups are left, the one closer in retention time to the expected value is chosen.

Two of the patterns extracted for $[PC \ 32:2 + H]^+$ are shown in Figure 1: the first (to the left) is relative to an unlabelled sample, whereas the second one (to the right) is relative to a labelled sample (99% 13 C labelling).

In this simple case, the difference is straightforward: in the labelled sample the most intense signal is shifted 40 mass units upwards with respect to the monoisotopic peak, indicating that the most abundant species is the one where all 40 carbon atoms have been replaced by the labelling isotope, ¹³C.

4.3 Isotopic Pattern Analysis

The extracted patterns are then analysed by the function find_abundance, which takes each of them and finds the best theoretical pattern that reproduces it by iteratively changing the relative abundance of the labelling isotope.

In the example, this information can be achieved through:

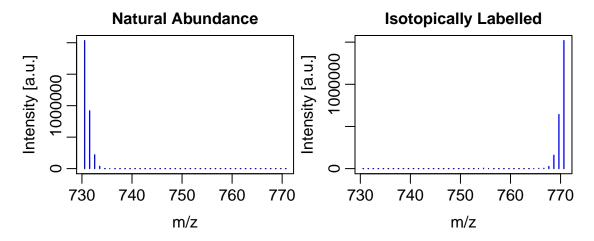


Figure 1: An example showing two of the patterns extracted from the experimental data; to the left is an unlabelled sample, to the right a labelled sample (99% ¹³C).

The output is a class *labelling* object, which is a list containing the results of the isotopic pattern analysis:

```
attributes(fitted_abundances)
```

The estimated percentage abundances of the labelling isotope are in "best_estimate" (numbers between 0 and 100), the standard errors from the fitting procedure are in "std_error", whereas the percentage deviations between best fitted and experimental patterns are in "dev_percent". The m/z values are in "x_scale", while in "y_exp" are their normalised intensities and in "y_theor" are those of the best fitted theoretical patterns. The differences between the two are in "residuals".

The find_abundance function performs the following:

- 1. It takes each experimental pattern and normalises it to its highest signal, set to 100;
- 2. If the initial estimates are not provided, it looks for the m/z position of the most intense signal, and uses it to get a first rough estimate of the X isotopic abundance;
- 3. A fitting procedure ensues, whereby the single unknown variable "isotopic abundance of X" is iteratively changed starting from its initial estimate, and the resulting theoretical patterns are compared to the experimental pattern. The final value is the one that minimizes the sum of squares of the difference between the two normalised patterns. In order to account for noise, a weight is given to the signals, proportional to the square root of their intensities. The theoretical patterns are computed using the ecipex R package, which exploits Fourier transforms of simplex-based elemental models.

5 Overview of the Results

There are a number of ways to look at and save the results of the isotopic pattern analysis:

1. The generic function summary allows to quickly glance at the estimated percentage abundances:

```
summary(fitted_abundances)
```

```
##
                      C12_Sample_1 C12_Sample_2 C12_Sample_3 C12_Sample_4
## Best Estimate [%]
                       1.074287762 1.080694206
                                                1.080582514
                                                             1.080180140
                       0.006673445 0.003860987
                                                0.004485302
## Standard Error [%]
                                                              0.005093723
##
                      C13_Sample_1 C13_Sample_2 C13_Sample_3 C13_Sample_4
## Best Estimate [%]
                      98.942141509 98.947395995 98.936148771 98.945395245
## Standard Error [%]
                      0.008717031 0.008215919 0.006310484
                                                             0.007579631
```

In this example, the average value is 1.08% for unlabelled samples (close to the natural 13 C abundance, 1.07%), whereas it is 98.94% for labelled samples.

2. The generic function plot can be used to produce three types of plots, depending on the parameter "type". By default (type="patterns") a series of plots is returned, one for each sample, showing the normalised experimental patterns superimposed to their fitted theoretical patterns:

```
plot(fitted_abundances, type="patterns", saveplots=F)
```

Two of the plots obtained in this example are in Figure 2.

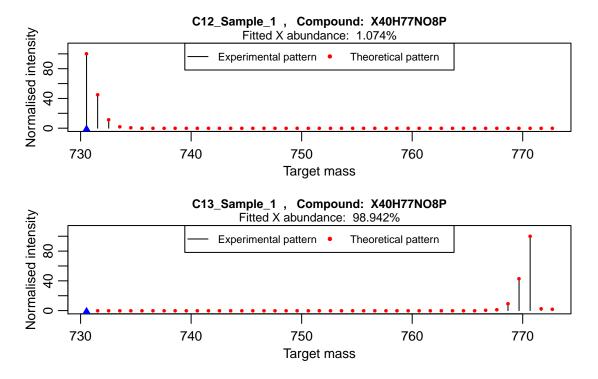


Figure 2: Graphical summary of the isotopic pattern analysis for an unlabelled (top) and a labelled (bottom) sample.

If "type" is set to "residuals", the residuals are plotted:

plot(fitted_abundances, type="residuals", saveplots=F)

This is shown in Figure 3.

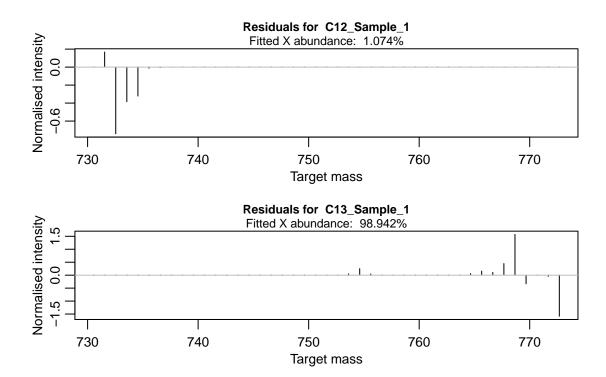


Figure 3: Plot of the residuals for an unlabelled (top) and a labelled (bottom) sample.

Finally, with type="summary", a summary plot with the estimated percentage abundances is provided (see Figure 4).

If the parameter "saveplots" is set to TRUE, the plots are saved as a *.pdf file in the working directory.

3. The save_labelling function allows to save the results to a *.csv file in the working directory:

```
save_labelling(fitted_abundances)
```

For each sample, this file reports:

- a. The estimated percentage abundance of the labelling isotope;
- b. The related standard error;
- c. The percentage deviation between theoretical and experimental isotopic patterns;
- d. The outcome message from the fitting procedure.

6 Conclusion

In this document we introduced the main functions of the *IsotopicLabelling* R package, explaining their basic working principles. Using the provided data set, we illustrated how to use the package in practice and how to quickly assess the results, with the hope that this package will prove to be a useful tool to researchers dealing with labelling experiments.

After having discussed the principles behind this package, our final aim here is to provide a compact summary that may be used as a script for analysing LC-MS data relative to labelling experiments.

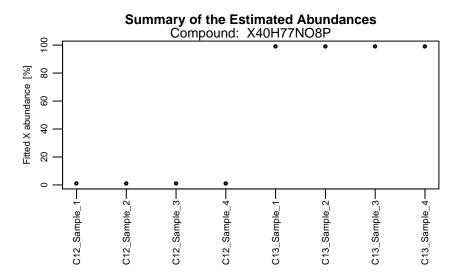


Figure 4: Graphical summary of the estimated percentage abundances and related standard errors, following the isotopic pattern analysis.

```
# Load the package
library("IsotopicLabelling")
# Load the xcmsSet object
data(xcms_obj)
# Convert the object into the required data frame
peak_table <- table_xcms(xcms_obj)</pre>
# Process the data
fitted_abundances <- main_labelling(peak_table, compound="X40H77N08P", labelling="C",
                                   mass_shift=0.05, RT=285, RT_shift=20,
                                    chrom_width=7, initial_abundance=NA)
# Quickly look at the results
summary(fitted_abundances)
# Plot the patterns
plot(fitted_abundances, type="patterns", saveplots=F)
# Plot the residuals
plot(fitted_abundances, type="residuals", saveplots=F)
# Plot the overall results
plot(fitted_abundances, type="summary", saveplots=F)
# Save the results to a *.csv file
save_labelling(fitted_abundances)
```